



# PRMEDICAL

PROMOTING INNOVATION

# Detecto

## Distal cap

## Instructions for use

### INTENDED USE

Detecto is intended for attachment to the distal end of the endoscope in order to facilitate endoscopic therapy, specifically improving visibility during colonoscopy screening procedures. Detecto can also be used to maintain adequate depth of the visual field of the endoscopy, gently stretch the internal folds of the colon, facilitate the advancement of the instrument through the gastrointestinal tract, stabilise the tip of the instrument inside the colon to offer more stable vision.

### ASSEMBLY

- 1 - Insert the distal end of the colonoscope into the lower cavity of Detecto
- 2 - Fully insert the tip of the endoscope until it touches the ring
- 3 - Make sure that the device is properly assembled. If Detecto is properly assembled, 2mm of spacer cap will protrude from the tip of the colonoscope, as shown in the figure.

It is possible to wet Detecto to facilitate the insertion of larger diameter colonoscopes into the device.

### CONTRAINDICATIONS

Do not use in patients known to have acute ulcerative colitis or colonic stenosis. The contraindications include those specific to any endoscopic procedure.

### PRECAUTIONS

Read the instructions carefully.

Before fitting Detecto onto the tip of the colonoscope, make sure that the end of the instrument and the device are in perfect condition and show no signs of damage.

Check that the outer diameter of the endoscope is compatible with the inner diameter of Detecto. See the following table for details of compatibility:

Commercial name	Ref	Box	Endoscope diameter
Detecto	PRM01	50pcs	12.5 - 13.5 mm



When fitting the Detecto device onto the colonoscope, hold the colonoscope as close as possible to the distal end in order to minimise the force applied to the flexible part of the colonoscope.

Do not apply lubricant gel to the colonoscope before attaching Detecto as this could cause the device to become detached from the colonoscope during use.

Detecto can be used on patients who have been prescribed a colonoscope, particularly members of the 50-74 age-group undergoing screening.

Detecto is a non-sterile, disposable class I device. Do not attempt to recondition or resterilise the device. PRMEDICAL did not design the device, nor is it intended, for reuse. The performance of these operations on this disposable medical device poses a risk to the safety of patients (compromising the integrity of the device, cross contamination, infection, for example).

The device is latex free. Keep the safety information label. In the event of serious incident relating to the device, please report this to PRMEDICAL and to the competent authority.

### DESCRIPTION OF SYMBOLS



See the instructions



Catalogue reference number



Batch code



Date of manufacture (mm) and (yy)



Disposable product



Non-sterile product



CE mark



Product Latex free



Manufacturer PRmedical s.r.l.

Via Castellana 16i, 30174, Venice Italy  
+39 351 9206411 info@prmedical.it




Medical Device



UDI code



Keep dry

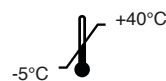

-5°C +40°C  
Temperature restrictions for storage and transport



5% 80%  
Humidity restrictions for storage and transport



Keep away from sunlight



### PRODUCT DISPOSAL

The product poses a potential biological risk of cross contamination after use. Handle the device and dispose of it in compliance with approved medical practice and the local, regional and national laws and regulations in force.

MADE IN ITALY